



HARTMANN Inco System



HARTMANN Inco System



The HARTMANN Inco System – in four colours.

The four colours of the HARTMANN incontinence products enable you to identify the absorption capacity of the individual products at a glance. The colour of the foil backing classifies the product. You know immediately what the product can do.

The four colours stand for the four degrees of incontinence. Yellow is for slight, green for moderate, blue for severe and violet for very severe incontinence. It's that simple. The colour coding appears not only on the product itself, but also on the inner and outer packaging.



Yellow for slight incontinence.



Green for moderate incontinence.



Blue for severe incontinence.



Violet for very severe incontinence.

The colour coding used for the HARTMANN Inco System brings you a number of benefits. It is easier to select the right product because you can immediately identify the product's performance from the colour of the foil backing. The practical information about product category and absorbency, shown on the packaging, enables quick and organised storage.

MoliNea® couche pads



Rectangular pads, with absorbent core of cellulose fluff.

Product description

An absorbent rectangular shaped pad for discreet protection against slight incontinence. The pad has a long adhesive strip which keeps the product reliably and firmly in place.

Uses

These products are designed for users with slight incontinence. Ideally to be used with Net Pants or snug fitted underwear.

Presentations



11 x 35 cm

MoliNea couche pads, size 2

No. 148520/6
224 (4 x P 56)



11 x 35 cm

MoliNea couche pads, size 3

No. 148530/6
224 (4 x P 56)



MoliMed®

Premium



In a compact design MoliMed is particularly thin and inconspicuous when worn.

Product description

A body shaped absorbent pad for discreet protection for slight incontinence. The pad has a long adhesive strip which keeps the product reliably and firmly in place. There are four absorbency levels available.

Uses

These products are designed for users with slight incontinence e.g. stress incontinence. Ideally to be used with Net Pants or snug fitted underwear.

Presentations



27 x 10 cm

MoliMed Premium micro

No. 168626/7
336 (24 x P 14)



27 x 14 cm

MoliMed Premium mini

No. 168636/7
252 (18 x P 14)



34 x 17 cm

MoliMed Premium midi

No. 168646/7
168 (12 x P 14)



45 x 21 cm

MoliMed Premium maxi

No. 168656/7
168 (12 x P 14)



MoliMed® for men



MoliMed for men protect is especially designed to cover the whole urogenital area and thus provides secure allround protection.

Product description

This is a specifically designed pad to fit the male anatomy. It offers discreet protection for slight incontinence. A wide adhesive strip on both versions keeps the product reliably and firmly in place.

Uses

These products are designed for male patients with slight incontinence e.g. after prostate surgery.

Available in two different versions:

MoliMed active – a pouch shaped pad for ultimate discreetness

MoliMed protect – a V shaped pad with inner anti-leak cuffs to cover the urogenital area for all round protection.

Presentations



13 x 15 cm

MoliMed for men active

No. 168600/7
168 (12 x P 14)



34 x 18 cm

MoliMed for men protect

No. 168707/7
168 (12 x P 14)



MoliForm®

Premium



MoliForm is available in four absorption capacities which makes it an ideal system for the management of all degrees of incontinence.

Product description

These are anatomically shaped pads with absorption capacities from light to very severe. The pads have integrated wetness indicators to show when the products needs changing.

Uses

MoliForm is available in four absorption capacities which makes it an ideal system for the management of all levels of incontinence. These products are chosen by the degree of absorbency required not size. MoliForm pads should be kept in place with appropriately sized Net Pants.

MoliPants®

soft



Elasticated mesh-pants for holding incontinence pads in place.

Product description

Soft textile fixation pants, with fine meshed material which is especially kind to the skin, with extra long leg for a neater fit.

Uses

The MoliPants soft are recommended for the perfect and secure fit of the shaped pad. The user's waist/hip should be measured as MoliPants soft are selected according to size. Size identification details are present on the inside of the waistband and also by colour code. MoliPants soft are re-useable, they can be washed up to 50 times at 60 °C and tumble dried at a low heat. Washing instructions can be found on inside of the waistband.

Presentations



62 x 32,5 cm

MoliForm Premium normal

No. 168022/6
112 (4 x P 28)



62 x 32,5 cm

MoliForm Premium plus

No. 168222/6
112 (4 x P 28)



69 x 36 cm

MoliForm Premium extra

No. 168322/6
112 (4 x P 28)



69 x 36 cm

MoliForm Premium super plus

No. 168922/6
112 (4 x P 28)



Presentations



yellow



Small



45 – 80 cm

MoliPants soft Small

No. 947790/2
25 (1 x P 25)



blue



Medium



60 – 100 cm

MoliPants soft Medium

No. 947791/2
25 (1 x P 25)



brown



Large



80 – 120 cm

MoliPants soft Large

No. 947792/2
25 (1 x P 25)



green



Large



100 – 160 cm

MoliPants soft X-Large

No. 947793/2
25 (1 x P 25)



grey



Large



140 – 180 cm

MoliPants soft XX-Large

No. 947794/2
25 (1 x P 25)



MoliPants® Economy



Elasticated net pants for holding incontinence pads in place.

Product description

Elasticated net fixation pants.

Uses

The MoliPants Economy are recommended for a more secure fit of shaped pads. The user's waist/hip should be measured as MoliPants Economy are selected according to size. Size identification is by the colour code on the waistband. These fixation pants are designed to be disposable.

Presentations



blue

60 – 100 cm

MoliPants Economy Medium

No. 947730/2
50 (1 x P 50)



brown

80 – 120 cm

MoliPants Economy Large

No. 947731/2
50 (1 x P 50)



green

100 – 180 cm

MoliPants Economy X-Large

No. 947738/1
50 (1 x P 50)

MoliCare® Premium Air Active, extra



The breathable non-woven side panels, covering the hip area makes the all-in-one incontinence briefs MoliCare Premium Air Active extra pleasant to wear.

Product description

This is an All In One incontinence brief (taped pad) with non-woven side panels enabling the skin to breathe. This pad has an integrated wetness indicator which shows when the product needs changing and elasticated adhesive tapes to allow for a greater flexibility.

Uses

Ideal for more heavily dependent patients, those at high risk of pressure sores, confused/agitated users or those with intractable diarrhoea/bowel management problems. The user's waist/hip should be measured as MoliCare Premium extra is selected according to size. The tapes should ideally be fixed onto the non woven side panels to ensure correct fitting of the product.

Presentations



SZ XS

40 – 60 cm

MoliCare Premium extra, Extra Small

No. 169960/5
112 (4 x P 28)



SZ 1

60 – 90 cm

MoliCare Premium extra, Small

No. 169961/5
84 (3 x P 28)



SZ 2

90 – 120 cm

MoliCare Premium extra, Medium

No. 169962/5
84 (3 x P 28)



SZ 3

120 – 150 cm

MoliCare Premium extra, Large

No. 169963/5
84 (3 x P 28)



SZ 4

150 – 175 cm

MoliCare Premium extra, Extra Large

No. 169964/5
56 (4 x P 14)



MoliCare® Classic extra



All-in-one incontinence briefs, with Dry Plus fluid distribution layer and wetness indicator.

Product description

This is an All In One incontinence brief (taped pad) with foil backed side panels and adhesive tapes. These pads have an integrated wetness indicator which shows when the product needs changing.

Uses

Ideal for more heavily dependent patients, those at high risk of pressure sores, confused/agitated users or those with intractable diarrhoea/bowel management problems.

The user's waist/hip should be measured as MoliCare Classic is selected according to size. The tapes should ideally be fixed on the foil side panels to ensure correct fitting of the product.

Presentations



SZ 1 50 – 80 cm

MoliCare Classic extra, Small
No. 169347/6
84 (3 x P 28)



SZ 2 70 – 120 cm

MoliCare Classic extra, Medium
No. 169547/6
84 (3 x P 28)



SZ 3 100 – 150 cm

MoliCare Classic extra, Large
No. 169747/6
84 (3 x P 28)



SZ 4 140 – 175 cm

MoliCare Classic extra, X-Large
No. 169944/6
56 (4 x P 14)



MoliCare® Premium Air Active super



The breathable non-woven side panels, covering the hip area makes the all-in-one incontinence briefs MoliCare Premium Air Active super pleasant to wear.

Product description

This is an All In One incontinence brief (taped pad) with non-woven side panels enabling the skin to breathe. This pad has an integrated wetness indicator which shows when the product needs changing and elasticated adhesive tapes to allow for a greater flexibility.

Uses

Ideal for patients who require a higher absorbency pad i.e. clients who have very severe incontinence.

It is also ideal for overnight use. Recommended for patients who are dependent, confused/agitated users or those with intractable diarrhoea/bowel management problems.

The user's waist/hip should be measured as MoliCare Premium super is selected according to size. The tapes should ideally be fixed onto the non woven side panels to ensure correct fitting of the product.

Presentations



SZ 1 60 – 90 cm

MoliCare Premium Air Active super, Small
No. 169971/5
84 (3 x P 28)



SZ 2 90 – 120 cm

MoliCare Premium Air Active super, Medium
No. 169972/5
84 (3 x P 28)



SZ 3 120 – 150 cm

MoliCare Premium Air Active super, Large
No. 169973/5
84 (3 x P 28)



MoliCare® Super Plus



Due to the extra large absorbent pad, the all-in-one incontinence briefs MoliCare Super Plus features a high fluid retention capacity.

Product description

This is an All In One incontinence brief (taped pad) with foil backed side panels and adhesive tapes. This pad has an integrated wetness indicator which shows when the product needs changing.

Uses

Ideal for patients who require a higher absorbency i.e. clients who have very severe incontinence. This pad has a high absorption rate which is of particular benefit when caring patients who experience high urinary volume outputs. It is also ideal for extended overnight use. Recommended for patients who are dependent, confused/agitated users or those with intractable diarrhoea/bowel management problems. The user's waist/hip should be measured as MoliCare super plus is selected according to size. The tapes should ideally be fixed on the foil side panels to ensure correct fitting of the product.

Presentations



SZ 1 50 – 80 cm

MoliCare Super Plus Small
No. 169470/5
56 (4 x P 14)



SZ 2 70 – 120 cm

MoliCare Super Plus Medium
No. 169670/5
56 (4 x P 14)



SZ 3 100 – 150 cm

MoliCare Super Plus Large
No. 169870/5
56 (4 x P 14)



MoliFlex®

Premium



Thanks to the belt fastening MoliFlex Premium incontinence pads can be securely held in place for maximum comfort.

Product description

This is a textile backed, belted incontinence pad. The belt fastening can be adjusted and readjusted thanks to the resealable velcro tapes. This pad also has an integrated wetness indicator which shows when the product needs changing.

Uses

Ideal for more heavily dependent patients, those at high risk of pressure sores, confused/agitated users or those with intractable diarrhoea/bowel management problems. MoliFlex is available in three absorption capacities which makes it an ideal system for the management of all levels of incontinence. The user's waist/hip should be measured as MoliFlex is selected according to size.

Presentations



SZ 1 60 – 85 cm

MoliFlex Premium Small
No. 169110/0
60 (3 x P 20)



SZ 2 80 – 105 cm

MoliFlex Premium Medium
No. 169240/0
63 (3 x P 21)



SZ 3 100 – 135 cm

MoliFlex Premium Large
No. 169340/0
63 (3 x P 21)



SZ 2 80 – 105 cm

MoliFlex maxi Medium
No. 169250/0
63 (3 x P 21)



SZ 3 100 – 135 cm

MoliFlex maxi Large
No. 169350/0
63 (3 x P 21)



SZ 2 80 – 105 cm

MoliFlex Premium super Medium
No. 169260/0
63 (3 x P 21)



SZ 3 100 – 135 cm

MoliFlex Premium super Large
No. 169360/0
63 (3 x P 21)

MoliCare® Mobile



MoliCare Mobile are reliable absorbent pants. They are inconspicuous when worn and can be pulled on and off like normal underwear.

Product description

This is a textile backed pull-up pants. This pad has an integrated wetness indicator which shows when the product needs changing, is discreet and can be fitted and removed like normal underwear.

Uses

Ideal for patients with moderate incontinence, but who are mobile and non-dependent.

This product is useful when toilet training and for those patients who are experiencing memory loss e.g. in the early stages of Alzheimer's.

The user's waist/hip should be measured as MoliCare Mobile is selected according to size.

Presentations



XS
Small



SZ XS

45 – 70 cm

MoliCare Mobile Extra Small

No. 915840/0
56 (4 x P 14)



S
Small



SZ 1

60 – 90 cm

MoliCare Mobile Small

No. 915841/0
56 (4 x P 14)



M
Medium



SZ 2

80 – 120 cm

MoliCare Mobile Medium

No. 915842/0
56 (4 x P 14)



L
Large



SZ 3

100 – 150 cm

MoliCare Mobile Large

No. 915843/0
56 (4 x P 14)



XL
Large



SZ 4

130 – 170 cm

MoliCare Mobile Extra Large

No. 915844/0
56 (4 x P 14)



MoliNea® L



Underpads with absorbent core of cellulose layers.

Product description

A waterproof under pad with an absorbent core, consisting of layers of cellulose.

Uses

Suitable for use on its own, or as additional bed protection for medical procedures.

Presentations



light
8 layers

57 x 60 cm

MoliNea L

No. 159926/6
100 (100 x P 1)



light
8 layers

57 x 90 cm

MoliNea L

No. 159933/6
100 (100 x P 1)



MoliNea® plus



Underpads with absorbent core of cellulose fluff.

Product description

A waterproof underpad with an absorbent core of cellulose fluff with a stay dry non-woven top layer.

Uses

Suitable for use on its own, or as additional bed protection for medical procedures.

Presentations



high absorbency

60 x 60 cm

MoliNea plus
No. 161300/6
100 (100 x P 1)



high absorbency

60 x 75 cm

MoliNea plus
No. 161540/6
100 (100 x P 1)



high absorbency

60 x 90 cm

MoliNea plus
No. 161600/6
100 (100 x P 1)

MoliNea® textile



Reusable Underpads.

Product description

A waterproof yet breathable underpad with a textured non-slip underside. The blue side is placed onto the mattress and the aqua side to the users skin. The underpad can be re-washed up to 300 times at 100 °C and is tumble dryable.

Uses

Suitable for use on its own, or as an additional bed protection.

Presentations



75 x 85 cm

MoliNea textile
No. 155801/2
30 (3 x P 10)



85 x 90 cm

MoliNea textile
No. 155802/2
30 (3 x P 10)



75 x 85 cm*

MoliNea textile
No. 155805/2
30 (3 x P 10)


* With tuckable flaps of 2 x 50 cm

All of the above items can be purchased as an individual item or in a pack of 10 depending on requirements.




MoliForm®


Fitting Guidelines – Standing

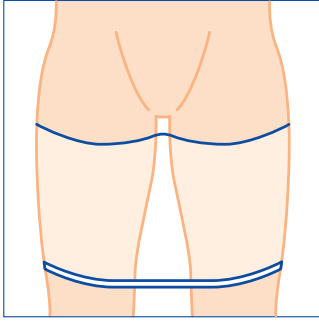


Step 1
Pull up MoliPants, roll waistband down to lower thigh

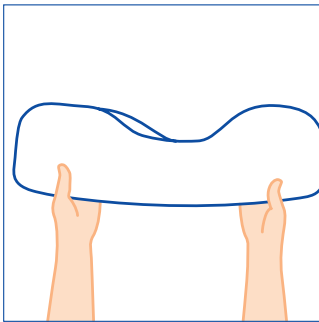


Step 2
Fold pad in half lengthways

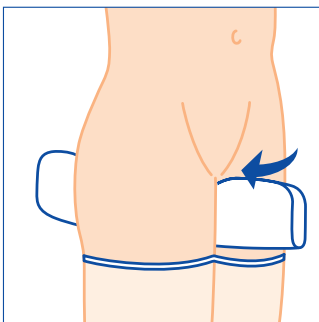

Step 3
Insert pad from front to back, and pull up tight into the groin

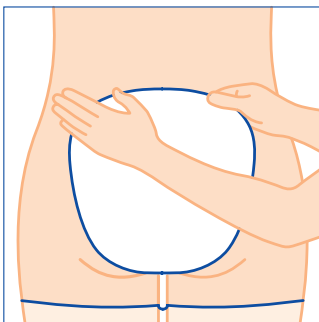


Step 4
Unfold and centre pad at the back, secure by pulling up MoliPants at rear

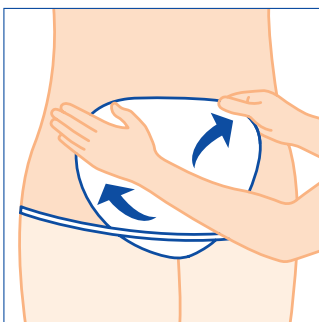


Step 5
Unfold and centre pad at the front, secure by pulling up MoliPants at front

MoliForm®


Fitting Guidelines – Lying Down

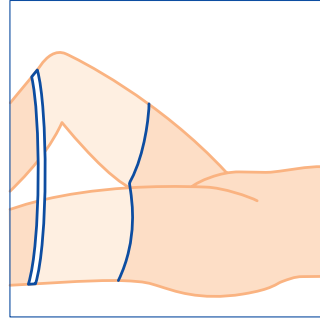


Step 1
Pull up MoliPants, roll waistband down to lower thigh and secure over knee

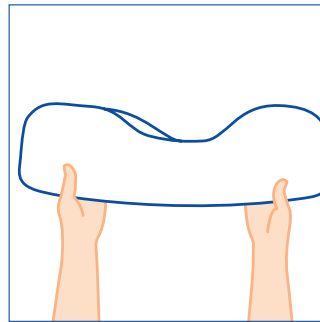


Step 2
Fold pad in half lengthways

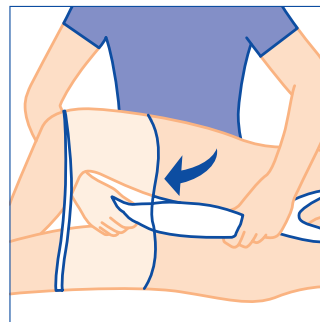

Step 3
Insert pad from front to back, and pull up tight into the groin

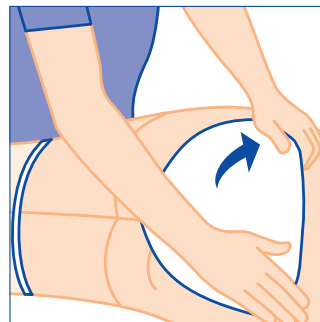


Step 4
Unfold and centre pad at the back, secure by pulling up MoliPants at rear

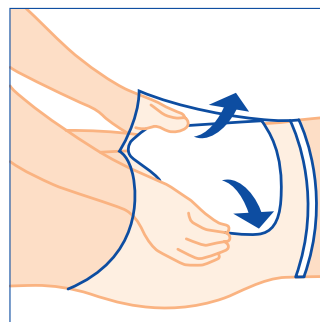


Step 5
Unfold and centre pad at the front, secure by pulling up MoliPants at front

MoliCare® extra


Fitting Guidelines – Standing

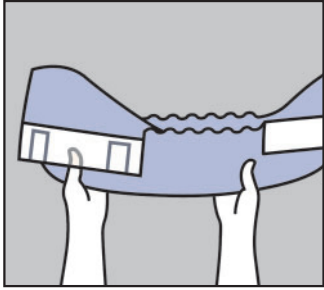


Step 1
Fold pad in half lengthways

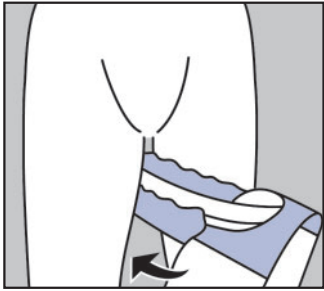

Step 2
Insert pad from front to back,
and pull up tight into groin

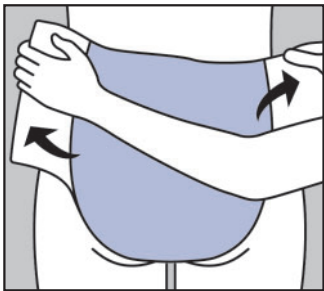

Step 3
Unfold and centre pad at the
back

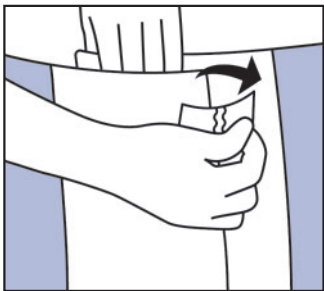

Step 4
Unfold and centre pad at the
front. Stretch and apply upper
tapes at downward angle




Step 5
Stretch and apply lower tapes at
upper angle



Step 6
Note: The tapes should sit over
the broadest part of hips, *not*
meeting in the middle

MoliCare® extra


Fitting Guidelines – Lying Down

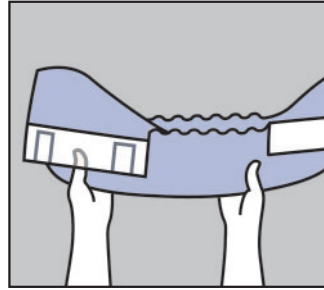

Step 1
Fold pad in half lengthways




Step 2
Insert pad from front to back,
and pull up tight into groin

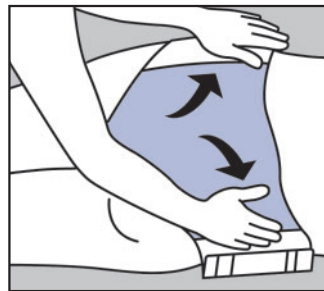


Step 3
Unfold and centre pad at the
back



Step 4
Unfold and centre pad at the
front. Stretch and apply upper
tapes at downward angle

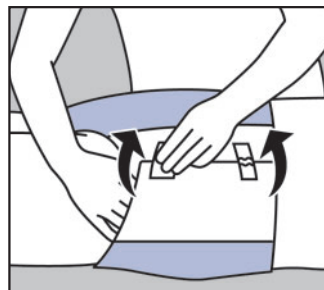

Step 5
Stretch and apply lower tapes at
upper angle

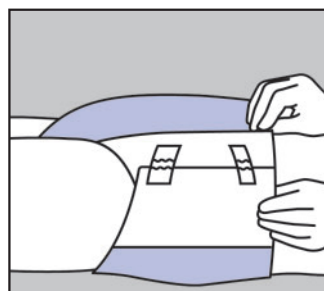

Step 6
Note: The tapes should sit over
the broadest part of hips, *not*
meeting in the middle



PAUL HARTMANN AG
Paul-Hartmann-Straße 12
89522 Heidenheim
P.O. Box 14 20
89504 Heidenheim
Phone +49-7321-36-0
Fax +49-7321-36-3636
www.hartmann.info

Germany



Corium
Grattan House
City Junction Business Park
Malahide Road
Dublin 17
Telephone: +353 1 866 1200
Facsimile: +353 1 877 4199
Web: www.corium.ie
NI Customers: 0845 601 7214

Ireland